

1

USAARL Report No. 93-23

AD-A269 967



Laboratory Test and Evaluation of the Essex Model 60C-0037-2 Molecular Sieve Oxygen Generating System

By

Brendan E. Squire
James E. Bruckart
Martin D. Quattlebaum
Leslie W. Johnston



Impact, Tolerance, and Protection Division

June 1993

93-22606



Approved for public release; distribution unlimited.

United States Army Aeromedical Research Laboratory
Fort Rucker, Alabama 36362-5292

Notice

Qualified requesters

Qualified requesters may obtain copies from the Defense Technical Information Center (DTIC), Cameron Station, Alexandria, Virginia 22314. Orders will be expedited if placed through the librarian or other person designated to request documents from DTIC.

Change of address

Organizations receiving reports from the U.S. Army Aeromedical Research Laboratory on automatic mailing lists should confirm correct address when corresponding about laboratory reports.

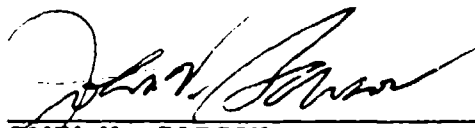
Disposition

Destroy this document when it is no longer needed. Do not return it to the originator.

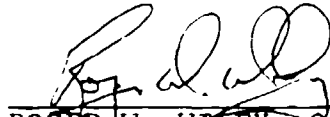
Disclaimer

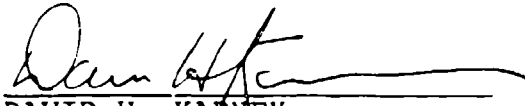
The views, opinions, and/or findings contained in this report are those of the author(s) and should not be construed as an official Department of the Army position, policy, or decision, unless so designated by other official documentation. Citation of trade names in this report does not constitute an official Department of the Army endorsement or approval of the use of such commercial items.

Reviewed:


JOHN V. BARSON
LTC, MC, SFS
Director, Biodynamics
Research Division

Released for publication:


ROGER W. WILEY, O.D., Ph.D.
Chairman, Scientific
Review Committee


DAVID H. KARNEY
Colonel, MC, SFS
Commanding

Unclassified

SECURITY CLASSIFICATION OF THIS PAGE

REPORT DOCUMENTATION PAGE				Form Approved OMB No. 0704-0188	
1a. REPORT SECURITY CLASSIFICATION Unclassified			1b. RESTRICTIVE MARKINGS		
2a. SECURITY CLASSIFICATION AUTHORITY			3. DISTRIBUTION/AVAILABILITY OF REPORT Approved for public release, distribution unlimited		
2b. DECLASSIFICATION/DOWNGRADING SCHEDULE					
4. PERFORMING ORGANIZATION REPORT NUMBER(S) USAARL Report No. 93-23			5. MONITORING ORGANIZATION REPORT NUMBER(S)		
6a. NAME OF PERFORMING ORGANIZATION U.S. Army Aeromedical Research Laboratory		6b. OFFICE SYMBOL (if applicable) SGRD-UAD-IE	7a. NAME OF MONITORING ORGANIZATION U.S. Army Medical Research and Development Command		
6c. ADDRESS (City, State, and ZIP Code) P.O. Box 620577 Fort Rucker, AL 36362-0577			7b. ADDRESS (City, State, and ZIP Code) Fort Detrick Frederick, MD 21702-5012		
8a. NAME OF FUNDING/SPONSORING ORGANIZATION		8b. OFFICE SYMBOL (if applicable)	9. PROCUREMENT INSTRUMENT IDENTIFICATION NUMBER		
8c. ADDRESS (City, State, and ZIP Code)			10. SOURCE OF FUNDING NUMBERS		
			PROGRAM ELEMENT NO 0603807A	PROJECT NO 3M463807D836	TASK NO LC
					WORK UNIT ACCESSION NO. 201
11. TITLE (Include Security Classification) Laboratory Test and Evaluation Report of the Essex Model 60C-0037-2 Molecular Sieve Oxygen Generating System					
12. PERSONAL AUTHOR(S) Brendan E. Squire, James E. Bruckart, Martin Quattlebaum, and Leslie W. Johnston					
13a. TYPE OF REPORT Final		13b. TIME COVERED FROM _____ TO _____		14. DATE OF REPORT (Year, Month, Day) 1993 June	
				15. PAGE COUNT 48	
16. SUPPLEMENTARY NOTATION					
17. COSATI CODES			18. SUBJECT TERMS (Continue on reverse if necessary and identify by block number)		
FIELD	GROUP	SUB-GROUP			
01	03				
14	02				
			Electromagnetic compatibility, test and evaluation, aeromedical equipment		
19. ABSTRACT (Continue on reverse if necessary and identify by block number)					
<p>The Essex Model 60C-0037-2 Molecular Sieve Oxygen Generating System was tested for environmental and electromagnetic interference/compatibility under the U.S. Army Program for Testing and Evaluation of Equipment for Aeromedical Operations. The tests were conducted using current military and industrial standards and procedures for environmental tests and electromagnetic interference/compatibility and human factors. The Essex Model 60C-0037-2 performed properly in the test environments. The oxygen concentration produced by the unit depends on the compressed air pressure, oxygen flow, and vent pressure. A separate oxygen concentration analyzer would be required to determine the oxygen concentration for patient use. The unit produces a significant amount of noise while operating.</p>					
20. DISTRIBUTION/AVAILABILITY OF ABSTRACT <input checked="" type="checkbox"/> UNCLASSIFIED/UNLIMITED <input type="checkbox"/> SAME AS RPT <input type="checkbox"/> DTIC USERS			21. ABSTRACT SECURITY CLASSIFICATION Unclassified		
22a. NAME OF RESPONSIBLE INDIVIDUAL Chief, Scientific Information Center			22b. TELEPHONE (include Area Code) 205-255-6907		22c. OFFICE SYMBOL SGRD-IAX-SI

Table of contents

SECTION		PAGE
1.	<u>EXECUTIVE DIGEST</u>	
1.1	Test objectives	1-1
1.2	Testing authority	1-2
1.3	Scope	1-2
1.4	Material description	1-2
1.5	Summary of findings	1-3
1.6	Conclusions	1-3
2.	<u>SUBTESTS</u>	
2.1	Initial inspection	2-1
2.2	Electrical safety evaluation	2-1
2.3	Pressure-flow performance	2-2
2.4	Human factors evaluation (laboratory)	2-2
2.5	Altitude (low pressure) test	2-3
2.6	Vibration test	2-4
2.7	High temperature test	2-5
2.8	Low temperature test	2-6
2.9	Humidity test	2-8
2.10	Electromagnetic characteristics test	2-9
2.11	Acoustical analysis	2-11

Accession For	
NTIS CRA&I	<input checked="" type="checkbox"/>
DTIC TAB	<input type="checkbox"/>
Unannounced	
Justification	
By	
Distribution	
Availability Codes	
Dist	Availability Codes Special
A-1	

3. SUPPORTING DOCUMENTATION

3.1	Photographic description	3-1
3.2	Detailed test information	3-2
3.3	Criteria, significant problems, and suggested improvements	3-26
3.4	References	3-28
3.5	Abbreviations	3-29
3.6	List of manufacturers	3-31

Section 1. Executive digest

There is a need to provide a continuous flow of oxygen for patient resuscitation on medical evacuation (MEDEVAC) aircraft. The present system of carrying bottled oxygen on MEDEVAC aircraft is logistically cumbersome, time critical, and potentially explosive in a combat environment. Molecular sieve oxygen generating systems (MSOGS) will eliminate supply, safety, and storage problems associated with using and securing pressurized oxygen bottles in the patient compartment of air ambulances.

The Army program for test and evaluation of aeromedical equipment uses existing military standards (MIL-STD) and collective professional expertise to test and evaluate selected medical equipment for use aboard Army aircraft. Equipment is tested to protect the patients, crew, and aircraft by eliminating problems due to environmental extremes, physical incompatibility with the aircraft, or electromagnetic interference with aircraft systems.

1.1 TEST OBJECTIVES

- 1.1.1 To determine if the MSOGS equipment is complete and operational per the manufacturer's operating instructions.
- 1.1.2 To ensure the electrical safety of the medical equipment.
- 1.1.3 To ensure the equipment will function as designed with various inlet air pressures and flow rates.
- 1.1.4 To ensure the safety of the operator, the patient, and the aircrew.
- 1.1.5 To assess design considerations which, potentially, could contribute to an operator error.
- 1.1.6 To determine if the MSOGS unit can function as designed in a low-pressure environment.
- 1.1.7 To determine the ability of the MSOGS equipment to withstand the vibrational stresses expected in a rotary-wing flight environment without degradation or malfunction.
- 1.1.8 To determine the ability of the MSOGS equipment to be stored and operated in a high-temperature environment.
- 1.1.9 To determine the ability of the MSOGS equipment to be stored and operated in a low-temperature environment.
- 1.1.10 To determine the ability of the MSOGS equipment to operate satisfactorily for short periods during exposure to highly humid conditions.

1.1.11 To assess the levels of electromagnetic emissions produced by the MSOGS equipment within selected frequency ranges.

1.1.12 To assess the minimum electromagnetic susceptibility levels of the MSOGS equipment within selected frequency ranges.

1.1.13 To assess the acoustical levels generated by the unit.

1.2 TESTING AUTHORITY

Research and Technology Work Unit Summary, dated 30 September 1991. Project number 3M463807D836, titled, Molecular Sieve Oxygen Generating Systems.

1.3 SCOPE

This test was conducted at the United States Army Aeromedical Research Laboratory (USAARL) at Fort Rucker, Alabama.

1.4 MATERIAL DESCRIPTION

A molecular sieve oxygen generator is an on-site oxygen generating machine. Coupled with a compressed air source and electrical power source, it separates oxygen from other gases. The separation is accomplished with an inert ceramic material that does not require replacement. The process is regenerative and should require little maintenance. Most oxygen generator systems use several ceramic filter beds to allow continuous operation. While one filter bed is "on-line" for the oxygen generation process, the other filter bed is cleaned and purged of nitrogen and other impurities. One is producing oxygen while another is purged of nitrogen.

The Essex* Model 60C-0037-2 Molecular Sieve Oxygen Generating System is designed to provide aviator breathing oxygen for crews in a high performance fixed wing aircraft. Oxygen can be delivered to the patient only when the generator is operating since there is no backup oxygen supply (BOS). The MSOGS requires a compressed air source and electrical power to produce oxygen. When installed in the UH-60 aircraft the engine bleed air will provide the compressed air and the aircraft electrical system will provide the electrical power to operate the unit. For the laboratory tests a commercial air compressor provided compressed air. Nitrogen is exhausted from the unit through a waste gas port on the unit.

* See list of manufacturers.

The Essex MSOGS filters compressed air and directs it to one of three Zeolite concentrator beds. The Zeolite material traps nitrogen molecules under pressure. Each bed is filled with compressed air and pressurized. Nitrogen molecules are trapped (filtered) in the Zeolite material and oxygen molecules are vented to an outlet port. When the bed is filled with nitrogen, the bed is depressurized and the nitrogen is purged. The oxygen concentration delivered by the unit depends on the compressed air inlet pressure, outlet pressure, oxygen flow rate, and altitude. A separate oxygen concentration analyzer is required to measure the concentration of oxygen produced by the unit.

1.5 SUMMARY OF FINDINGS

1.5.1 Electrical Safety Evaluation: All measurements were within acceptable limits. No unsafe qualities were found in the Essex MSOGS.

1.5.2 Pressure-flow Performance: The unit will function with a wide range of available inlet pressure and flow conditions. If sufficient compressed air pressure or flow is not available, the unit may produce air with less than 90 percent oxygen.

1.5.3 Human Factors Evaluation: The Essex MSOGS was found to be satisfactory in all applicable categories of the evaluation.

1.5.4 Environmental Tests: The Essex unit can be expected to perform in a variety of environmental conditions. Its performance was found to be satisfactory in all stages of the environmental testing. The requirements for environmental tests are established in MIL-STD-810E, Methods 500.2 (altitude), 514.3 (vibration), 501.2 (high temperature), 502.2 (low temperature), and 507.2 (humidity).

1.5.5 Electromagnetic Characteristics Tests: No unsafe qualities were found in the Essex MSOGS unit. The unit may be unsatisfactory in certain EMI sensitive environments. Some broadband radiated and conducted emissions exceeded the test limits.

1.5.6 Acoustical Analysis: The Essex MSOGS produces a significant amount of noise. Noise levels from 60 dB to 80 dB were measured from 1000 to 10,200 Hz.

1.6 CONCLUSIONS

The Essex Molecular Sieve Oxygen Generator performed properly in the test environments. The oxygen concentration produced by the unit depends on the compressed air pressure, oxygen flow, and vent pressure. A separate oxygen concentration analyzer would be required to determine the oxygen concentration for

patient use. The unit produces a significant amount of noise while operating.

Section 2. Subtests

2.1 INITIAL INSPECTION

2.1.1 Objective

To determine if the MSOGS unit is complete and operational for testing per the manufacturer's operating instructions.

2.1.2 Criteria

2.1.2.1 The physical inventory is conducted solely for investigation and documentation.

2.1.2.2 The MSOGS unit will be operated with the manufacturer's recommended power and compressed air provided. The unit will produce a gas of at least 90 percent oxygen concentration at the manufacturer's specified flow rate.

2.1.3 Test procedure

2.1.3.1 A complete physical inventory of the MSOGS unit will be completed per the manufacturer's equipment list.

2.1.3.2 An operational validation test of the MSOGS unit will be conducted per the manufacturer's operating instructions by USAARL's personnel.

2.1.4 Test findings

2.2 ELECTRICAL SAFETY EVALUATION

2.2.1 Objective

To ensure the electrical safety, by evaluation of case-to-ground resistance and case-to-ground current leakage, of the MSOGS unit.

2.2.2 Criterion

The MSOGS unit shall meet the standards established in NAF 99 for electrical safety of medical equipment.

2.2.3 Test procedure

The MSOGS unit shall be tested with a Neurodyne-Dempsey model 431F electrical safety analyzer* IAW the procedures described in Technical Bulletin (TB) Number 38-750-2. Checks for safety concerns such as case integrity, breaks in power cord insulation, and connectors shall be performed.

2.2.4 Test findings

An electrical safety evaluation of the control box for the unit could not be performed because the electrical supply connections were not supplied with the unit. A custom control box was fabricated at USAARL to allow testing of the unit. There were no safety concerns with the electrical connectors or power lines on the unit. Criterion partially met.

2.3 PRESSURE-FLOW PERFORMANCE (Laboratory)

2.3.1 Objective

To ensure the equipment will function as designed with inlet air flow rates from 50 to 200 percent of the flow rate specified by the manufacturer.

2.3.2 Criterion

The unit will produce oxygen of at least 90 percent concentration at the manufacturer's specified flow rate.

2.3.3 Test procedure

2.3.3.1 The MSOGS unit will be operated for at least 5 minutes with air flow rates from 50 to 200 percent (in 25 percent increments) of the value specified by the manufacturer.

2.3.3.2 The oxygen percentage and flow rate of the gas produced by the generator will be measured.

2.3.4 Test findings

The MSOGS unit consistently produced high concentration oxygen (greater than 90 percent) throughout a wide range of compressed air inlet pressure and output flow rates. Several test conditions are summarized:

<u>Inlet pressure (psi)</u>	<u>Oxygen flow (lpm)</u>	<u>Oxygen concentration</u>
25	15	90%
30	15	92.5%
40	15	94%

Criterion met.

2.4 HUMAN FACTORS EVALUATION (Laboratory)

2.4.1 Objectives

2.4.1.1 To assure the safety of the operator, the potential patient, and the aircrew.

2.4.1.2 To assess the design considerations which, potentially, could contribute to an operator error.

2.4.2 Criterion

The MSOGS unit must be rated satisfactory in all major categories of the evaluation. These include visual displays, controls, maintainability, conductors, fasteners, test points, test equipment, fuses and circuit breakers, labels and coding, and safety.

2.4.3 Test procedure

2.4.3.1 The evaluation will be conducted in a laboratory under fluorescent lighting and ambient room conditions.

2.4.3.2 The MSOGS unit will be operated according to prescribed instructions through its full range of functions.

2.4.4 Test findings

The unit was found to be satisfactory in all of the applicable criteria. Several criteria were not applicable due to the custom control box and installed characteristics of this device. Criterion met.

2.5 ALTITUDE (LOW PRESSURE) TEST [IAW MIL-STD-810E, METHOD 500.3]

2.5.1 Objective

To determine if the MSOGS unit can function as designed in a low-pressure environment.

2.5.2 Criterion

The MSOGS unit will produce oxygen of at least 90 percent concentration at the specified flow rate while exposed to an altitude equivalency of 15,000 feet above sea level.

2.5.3 Test procedure

2.5.3.1 A pretest performance check will be conducted to ensure proper operation of the MSOGS unit.

2.5.3.2 The altitude test was performed in a Tenney Engineering model 64S altitude chamber*. This test is based on MIL-STD-810E, Method 500.3. The MSOGS unit was operated on the floor of the chamber. Chamber pressure was decreased to 420 mmHg (15,000 ft equivalent altitude) over a 15-minute period, held constant for 60 minutes, then raised at 1500 fpm, to ambient conditions.

2.5.3.3 A posttest performance check will be conducted to ensure proper operation of the MSOGS unit after the exposure to low pressure.

2.5.4 Test findings

2.5.4.1 The Essex MSOGS met pretest performance criterion 2.2.2.1.

2.5.4.2 The unit operated properly and produced high concentration oxygen while at 15,000 feet equivalent altitude. Criterion met.

2.5.4.3 The unit operated properly after the altitude exposure. Criterion met.

2.6 VIBRATION TEST [IAW MIL-STD-810E, METHOD 514.4]

2.6.1 Objective

To determine the ability of the MSOGS unit to withstand the vibrational stresses expected in a rotary-wing environment without degradation or malfunction.

2.6.2 Criterion

The MSOGS unit will remain operational and be able to produce at least 90 percent oxygen at the specified flow rate while exposed to vibrational stresses.

2.6.3 Test procedure

2.6.3.1 A pretest performance check will be conducted to ensure proper operation of the MSOGS unit.

2.6.3.2 The vibration test was performed using an Unholtz-Dickey model TA115-40/CSTA vibration test system*. It is a single-axis system with an electromagnetic driver unit. The test consists of sinusoidal vibrations superimposed on random vibrations over a frequency range of 500 Hz, as shown below. These vibrations are derived from measurements taken on the floor under the copilot's seat in a UH-1 helicopter traveling at 120 knots. The reference spectrum breakpoints are from MIL-STD-810E, Method 514.4; reference spectrum levels are based on field measurements with a conservatism factor of 1.5. Independent tests were conducted in the X, Y, and Z-axes.

X, Y and Z-axis

duration: 60 minutes
broadband intensity: 2.116 G_{rms}
random vibration: initial slope : 99.00 dB/oct
5 Hz level: 0.00007 G_{sq}/Hz
100 Hz level: 0.0007 G_{sq}/Hz
300 Hz level: 0.0007 G_{sq}/Hz
500 Hz level: 0.00007 G_{sq}/Hz
final slope: -99.00 dB/oct

sinusoidal vibration:

0.9500 G_{pk} at 17.00 Hz
1.4000 G_{pk} at 34.00 Hz
0.9000 G_{pk} at 51.00 Hz
0.9000 G_{pk} at 68.01 Hz

The MSOGS unit was connected to the vibration table fixture, and its performance evaluated before, during, and after exposure to vibration.

2.6.3.3 A posttest performance check was conducted to ensure proper operation of the MSOGS unit.

2.6.4 Test findings

2.6.4.1 The pretest performance check met the criterion.

2.6.4.2 The vibration test met criterion 6.6.2.2.

2.6.4.3 The posttest performance check met the criterion
2.6.2.3.

2.7 HIGH TEMPERATURE TEST [IAW MIL-STD-810E, METHOD 501.3]

2.7.1 Objective

To determine the ability of the MSOGS unit to be stored and operated in a high-temperature environment.

2.7.2 Criteria

2.7.2.1 The MSOGS unit will generate at least 90 percent oxygen at the specified flow rate during the high-temperature operation check.

2.7.2.2 The MSOGS unit will generate at least 90 percent oxygen after the high-temperature storage cycle.

2.7.3 Test procedure

2.7.3.1 A pretest performance check was conducted to ensure proper operation of the MSOGS unit.

2.7.3.2 The high-temperature test was conducted in a Tenney Engineering model 2WUL-10107D walk-in controlled environment chamber*. This test is based on MIL-STD-810E, Method 501.3. For the high-temperature operation test, the MSOGS unit was operated in the environmental chamber. The chamber temperature was raised to 49°C and the humidity stabilized at a maximum of 20 percent RH within 15 minutes. The environmental control system is capable of regulating temperature within $\pm 2^\circ\text{C}$ and humidity within ± 5 percent RH. Temperature and humidity were held constant for 2 hours. At 30-minute intervals, the chamber door was opened briefly to check the oxygen percentage and flow rate. After the operational test, the MSOGS unit was allowed to return to ambient conditions over a 30-minute period.

2.7.3.3 A posttest performance check was conducted to ensure proper operation of the MSOGS unit.

2.7.3.4 The MSOGS unit was stored (not operated) at temperatures of 63°C for 1 hour, 71°C for 4 hours, then again at 63°C for 1 hour. The chamber and MSOGS unit were returned to ambient conditions over a 30-minute period.

2.7.3.5 A poststorage performance check was conducted to ensure proper performance of the MSOGS unit.

2.7.4 Test findings

2.7.4.1 The pretest performance check met criterion 2.7.2.1.

2.7.4.2 The Essex Model 60C-0037-2 performed normally during the high temperature test. Criterion met.

2.7.4.3 The posttest performance check met criterion 2.7.2.3.

2.7.4.4 The Essex Model 60C-0037-2 functioned properly after high temperature storage. Criterion met.

2.8 LOW TEMPERATURE TEST [IAW MIL-STD-810E, METHOD 502.3]

2.8.1 Objective

To determine the ability of the MSOGS unit to be operated and stored in a low-temperature environment.

2.8.2 Criteria

2.8.2.1 The MSOGS unit will produce at least 90 percent oxygen at the specified flow rate during the low-temperature operation check.

2.8.2.2 The MSOGS unit will produce at least 90 percent oxygen at the specified flow rate after the low-temperature storage cycle.

2.8.3 Test procedure

2.8.3.1 A pretest performance check was conducted to ensure proper operation of the MSOGS unit.

2.8.3.2 The MSOGS unit was placed on the floor of the environmental chamber and the temperature lowered to -25°C within 25 minutes. The environmental control system is capable of regulating temperature within 2°C. Humidity cannot be controlled in the chamber at freezing temperatures. The temperature was held constant for 2 hours. The chamber door was opened briefly every 30 minutes to check the percentage and flow rate of oxygen produced. The chamber temperature was raised to ambient temperature within a 30-minute period.

2.8.3.3 A posttest performance check was conducted to ensure proper operation of the MSOGS unit.

2.8.3.4 The MSOGS unit was "stored" in a nonoperational mode. The MSOGS unit was placed on the floor of the environmental test chamber and the temperature lowered to -46°C for 6 hours. The chamber was raised to ambient temperature over a 30-minute period.

2.8.3.5 A poststorage performance check was conducted to ensure proper operation of the MSOGS unit.

2.8.4 Test findings

2.8.4.1 The pretest performance check met criterion 2.8.2.1.

2.8.4.2 The Essex Model 60C-0037-2 performed normally during exposure to a cold environment. Criterion met.

2.8.4.3 The posttest performance check met criterion 2.8.2.3.

2.8.4.4 The MSOGS unit functioned properly after the low temperature storage test. Criterion met.

2.9 HUMIDITY TEST [IAW MIL-STD-810E, METHOD 507.3]

2.9.1 Objective

To determine the ability of the MSOGS unit to operate satisfactorily for short periods of time during exposure to highly humid conditions.

2.9.2 Criterion

The MSOGS unit will produce at least 90 percent oxygen at the specified flow while exposed to a high-humidity environment.

2.9.3 Test procedure

2.9.3.1 A pretest performance check was conducted to ensure the proper operation of the MSOGS unit.

2.9.3.2 The humidity test was conducted in a Tenney Engineering model ZWUL-10107D walk-in controlled environment chamber*. This test is based on MIL-STD-810E, Method 507.3. For the humidity test, the MSOGS unit was placed in operation on the floor of the environmental chamber. The chamber temperature was raised to a temperature of 30°C and a relative humidity of 95 percent within 25 minutes. Temperature and relative humidity were maintained for 4 hours. The environmental control system is capable of regulating temperature within $\pm 2^\circ\text{C}$ and humidity within ± 5 percent RH. At 45-minute intervals, the performance of the oxygen generator was checked. The chamber and the MSOGS unit were returned to ambient conditions before the posttest performance validation check was conducted.

2.9.3.3 A posttest performance check was conducted to ensure the proper operation of the MSOGS unit.

2.9.4 Test findings

2.9.4.1 The pretest performance check met criterion 2.9.2.1.

2.9.4.2 No failures were noted in performance checks during exposure to a high humidity environment. Criterion met.

2.9.4.3 The posttest performance check met criterion 2.9.2.3.

2.10 ELECTROMAGNETIC CHARACTERISTICS TEST [IAW MIL-STD-461C, Notice 2, and MIL-STD-462, Notice 5]

2.10.1 Objectives

2.10.1.1 To assess the maximum levels of radiated electromagnetic emissions produced by the MSOGS unit in the 14 kHz to 12.4 GHz frequency range.

2.10.1.2 To assess the tolerances of radiated electromagnetic susceptibility of the MSOGS unit within the 10 kHz to 10 GHz electric field.

2.10.2 Criteria

2.10.2.1 The MSOGS unit will not produce emissions in excess of the limits set forth in MIL-STD-461C, Notice 2.

2.10.2.2 The MSOGS unit will not malfunction when it is subjected to radiated emissions as specified in MIL-STD-461C, Notice 2.

2.10.2.3 The MSOGS unit shall not conduct emissions in excess of the limits set forth in MIL-STD-461C, Notice 2.

2.10.2.4 The MSOGS unit shall not malfunction when it is subjected to conducted emissions as specified in MIL-STD-461C, Notice 2.

2.10.3 Test procedure

2.10.3.1 The radiated emissions test was performed according to MIL-STD-462, Notice 5, Method RE02. The MSOGS unit was positioned on a wooden test stand inside the EMI chamber, 1 meter away from the receiving antennas. The antennas are mounted for both vertical and horizontal polarities and connected to EMI receivers. The MSOGS unit was connected to an air source through a tube extended outside the chamber. While the MSOGS unit was operating, the frequency spectrum (14 kHz to 12.4 GHz) was scanned for emissions.

2.10.3.2 The radiated susceptibility test was performed according to MIL-STD-462, Notice 5, Method RS03. The MSOGS unit was positioned on a wooden test stand inside the EMI chamber, 1 meter away from the transmitting antennas. The antennas are mounted for both vertical and horizontal polarities and connected to radio frequency (RF) transmitters. The MSOGS unit was connected to an air source through a tube extended outside the chamber. While the MSOGS unit was operating, it was monitored for faulty operation during exposures to fields of 20 V/m from 10 kHz to 10 GHz.

2.10.3.3 The conducted emissions tests will be performed according to MIL-STD-462, Notice 5, Methods CE01 and CE04. The MSOGS unit was placed on a grounded, copper-covered workbench. The top of the workbench is 1 meter from floor level, 1.37 meters long and 0.81 meters wide. Power was supplied via a pair of line impedance stabilization networks (LISN) and a test jig. The test jig is a wooden tray with two power receptacles and two slots to hold current probes in place around power supply conductors. While the MSOGS unit was operating, the frequency range (10 kHz to 50 MHz) was scanned for emissions conducted in the power cable from the MSOGS unit.

2.10.3.4 The conducted susceptibility spike test was performed according to MIL-STD-462, Notice 5, Method CS06, on a chemical resistant counter top. Power was supplied via a customized metal connection box. The connection box has two power receptacles and four banana jacks on its front panel. Connections to the individual power lines are made in series through the banana jacks. Transient spikes of 56 volts, 10 microseconds are generated with a Solar Electronics model 8282-1 transient pulse generator* and induced onto the power leads at the connection box banana jacks. The spikes are monitored with a Tektronix 2235 oscilloscope* connected to a power receptacle on the connection box. The MSOGS unit was plugged into the other receptacle on the connection box and placed in operation. It was observed for proper operation and visual displays while it was subjected to the power line spikes.

2.10.3.5 The conducted susceptibility test was performed according to MIL-STD-462, Notice 5, Method CS02. The MSOGS unit was placed on a grounded, copper-covered workbench. Radio frequency interference was induced on the power leads and measured at the MSOGS unit power cable. The frequency of the interference is incremented over the 50 kHz to 400 MHz range while the MSOGS unit was operated. It was observed for proper operation and visual displays while it was subjected to the radio interference on the power leads. Each frequency was held for 15 seconds.

2.10.4 Test findings

2.10.4.1 During the radiated emissions test a narrowband and broadband emissions exceeded the specification limits. The narrowband emission at 3.789 MHz exceeded the specification by 52.7 dB. Broadband emissions were detected from 0.014 - 625 MHz with 0.4 - 32.6 dB of failure. Criterion partially met.

2.10.4.2 The Essex 60C-0037-2 was not susceptible to radiated frequency interference. Criterion met.

2.10.4.3 Signal failures were detected during the conducted emissions tests. CE01 emissions were detected from 23.42 - 50

MHz with 0.1 - 7.9 dB of failure. CE04 emissions were detected from 0.325 - 45 MHz with 0.1 - 44.9 dB of failure. Criterion partially met.

2.10.4.4 The Essex 60C-0037-2 was not susceptible to radio frequency interference or test spikes during the conducted susceptibility tests. Criterion met.

2.11 ACOUSTICAL ANALYSIS

2.11.1 Objective

To determine the noise hazard associated with operating the Essex Model 60C-0037-2 Molecular Sieve Oxygen Generating System.

2.11.2 Criterion

The Essex MSOGS will not produce a noise hazard during normal operation.

2.11.3 Test procedure

The acoustical analysis was conducted in an open outdoor area at USAARL. A Nagra recorder* was positioned 20 ft. from the MSOGS unit and connected to two microphones positioned 4 feet from the unit at an angle 90 degrees away from each other. Acoustic samples were taken north, east, south, and west of the unit. The MSOGS unit exhaust port system was not constrained or muffled during testing.

2.11.4 Test findings

The Essex Model 60C-0037-2 produces a noise which is subjectively loud. The peak noise level was 82.9 dB. Noise levels from 60 dB to 80 dB were produced in the 1000 Hz to 10,200 Hz range. Additional testing when the unit is installed will be required to properly evaluate this criterion for operation in the aircraft. Criterion not evaluated.

Section 3. Supporting documentation

3.1 PHOTOGRAPHIC DESCRIPTION



3.2 DETAILED TEST INFORMATION

3.2.1 Electrical safety test

Electrical Safety Test Report Form

Nomenclature: Molecular Sieve Oxygen Generating System
Manufacturer: Essex
Model number: 60C-0037-2
Serial number: 0002
Military item number: None

Options installed: None

Date of test: May 92

Performance:

Grounding conductor resistance (milliohms): NA

Leakage current - Case to ground (microamperes):

unit off, grounded, normal polarity	NA
unit off, ungrounded, normal polarity	NA
unit off, ungrounded, reverse polarity	NA
unit on, grounded, normal polarity	NA
unit on, ungrounded, normal polarity	NA
unit on, ungrounded, reverse polarity	NA

MAXIMUM LIMITS:

ground resistance (milliohms):	150
current (microamperes)	
current (grounded, type A unit):	10
current (ungrounded, type A unit):	100
current (grounded, type B unit):	50
current (ungrounded, type B unit):	500

Comments on item setup or checks: A control box was not provided by the manufacturer and was fabricated at USAARL. Testing of this device is not applicable to a unit installed in the aircraft.

Comments on test run (including interruptions): No unsafe qualities were noted in the unit.

Comments on other data: None

3.2.2 Human factors evaluation

Human Factors Evaluation Report Form

Nomenclature: Molecular Sieve Oxygen Generating System
Manufacturer: Essex
Model number: 60C-0037-2
Serial number: 0002
Military item number: None

Options installed: None

Date of test: May 92

Item configuration during test: Item prepared for operation.

Checklist for HFE

RESULTS

VISUAL DISPLAYS:

NA

display type, format, content
location of displays
indicator lights
scalar displays
color coding
legends and labels
cathode ray tubes
counters
flags, go-no-go, center-null indicators

Comments: Visual displays not provided on the unit.

CONTROLS:

Satisfactory

location
characteristics of controls
labeling
control - display relationships

Comments:

TIME REQUIRED TO PREPARE FOR OPERATION (list in comment)

Comments:

MAINTAINABILITY:

Satisfactory

component location
component characteristics
rests and stands
covers, cases, access doors
handles
lubrication
component mounting
cord storage provisions
external accessibility
internal accessibility
list special tools required
list realistic inspection requirements
list realistic inspection intervals

Comments: Only item of user maintenance is replacement of air filter.

CONDUCTORS:

Satisfactory

binding and securing
length
protection
routing
conductor coding
fabrication
connectors

Comments: Control box not provided with the unit.

FASTENERS:

Satisfactory

access through inspection panel covers
enclosure fasteners
device mounting bolts and fasteners

Comments: None

TEST POINTS:

NA

general
location and mounting
test point labeling and coding

Comments: None.

TEST EQUIPMENT:

NA

general
equipment self-test
indicators (list in comments)
controls
positive indication of proper operation

Comments:

FUSES AND CIRCUIT BREAKERS:

Satisfactory

external accessibility
easy replacement or reset by operator

Comments: Manufactured control box was equipped with fuses.

LABELS AND CODING:

Satisfactory

placed above controls and displays
near or on the items they identify
not obscured by other equipment components
describe the function of the items they identify
readable from normal operating distance
conspicuous placards adjacent to hazardous items

Comments: None

SAFETY:

Satisfactory

manual
materials
fire and explosive protection
operator protection from mechanical hazards
patient protection from mechanical hazards
electrical safety (operator and patient)

Comments: No manuals supplied with the unit.

3.2.3 Altitude Tests

Altitude Test Report Form

Nomenclature: Molecular Sieve Oxygen Generating System
Manufacturer: Essex
Model number: 60C-0037-2
Serial number: 0002
Military item number: None

Options Installed: None

Date of test: 17 Aug 92

Item configuration during test:

Operating on the chamber floor and connected to the USAARL oxygen test apparatus (OTA) via portals in the chamber wall.

Performance test criteria:

The unit's ability to provide a continuous flow (16 lpm) of 90% (or higher) oxygen with an input pressure of 60 psi.

Ambient conditions outside chamber:

Temperature	77°F
Humidity	66% RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria)
all OK: Pass

Installation of item in test facility:

list connections to power	28 Vdc
list connections to simulators	Oxygen tester
list connections to dummy loads	None
list unconnected terminals	None

IN-test DATA

Time of test start: 1945

POSTTEST DATA

Posttest performance check:

(complete check of item and accessories)

Time of test end: 1145

Item functional (based on performance test
criteria): all OK: Pass

Deviation from pretest: None

Comments on item setup or checks:

Comments on test run (including interruptions): None

Comments on other data:

3.2.4 Vibration tests

Vibration test Report Form

Nomenclature: Molecular Sieve Oxygen Generating System
Manufacturer: Essex
Model number: 60C-0037-2
Serial number: 0002
Military item number: None

Options Installed: None

Date of test: 19 Aug 92

Item configuration during test: Attached to a custom
aluminum vibration table fixture.

Performance test criteria:

The unit's ability to provide a continuous flow (16
lpm) of 90% (or higher) oxygen with an input pressure
of 60 psi.

PRETEST DATA

Pretest performance check: All OK: Pass
Item functional (based on performance test criteria)
All OK: Pass

Installation of item in test facility:

list connections to power	28 Vdc
list connections to simulators	Oxygen test apparatus
list connections to dummy loads	None
list unconnected terminals	None

Ambient conditions

Temperature 77°F
Humidity 65% RH
Barometric pressure 1 atm

IN-TEST DATA

Data and performance checks during test:

Times and dates of test start:

X:9/19/92 1515 Y:9/18/92 1245 Z:9/20/92 1350

Time at first check:

X:1525 Y:1255 Z:1400

Item functional (based on performance test criteria):

All OK: Pass

Deviation from pretest: None

Time at second check:

X:1505 Y:1235 Z:1340

Item functional (based on performance test criteria):

All OK: Pass

Deviation from pretest: None

POSTTEST DATA

Time at test end:

X:1615 Y:1345 Z:1450

Posttest performance check: All OK: Pass

(complete check of item and accessories)

Item functional (based on performance test criteria):

All OK: Pass

Item intact: Yes

Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: None

3.2.5 High Temperature Operation test

High Temperature test (Equipment Operating) Report Form

Nomenclature: Molecular Sieve Oxygen Generating System
Manufacturer: Essex
Model number: 60C-0037-2
Serial number: 0002
Military item number: None

Options installed: None

Date of test: 19 Jun 92

Item configuration during test:

Operating on a wire test stand and connected to the
USAARL oxygen test apparatus via portals in the chamber
wall.

Performance test criteria:

The unit's ability to provide a continuous flow (16
lpm) of 90% (or higher) oxygen with an input pressure
of 60 psi.

Ambient conditions outside chamber:

Temperature	26°C
Humidity	79% RH
Barometric pressure	752.5 mmHg

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria):
All OK: Pass

Installation of item in test facility:

list connections to power	28 Vdc
list connections to simulators	OTA
list connections to dummy loads	none
list unconnected terminals	none
distance from north wall (meters)	0.56
distance from south wall (meters)	1.02
distance from east wall (meters)	1.57
distance from west wall (meters)	1.45
distance from ceiling (meters)	1.19
distance from floor (meters)	0.97

Time of test start: 1315

Performance checks during test:

First check:

Time: 1345
Temperature: $49^{\circ}\text{C} \pm 1^{\circ}\text{C}$
Humidity: $15\% \text{ RH} \pm 1\% \text{ RH}$
Barometric pressure: 752.5 mmHg
Item functional (based on performance criteria):
All OK: Pass
Deviation from pretest: None

Second check:

Time: 1415
Temperature: $49^{\circ}\text{C} \pm 1^{\circ}\text{C}$
Humidity: $15\% \text{ RH} \pm 1\% \text{ RH}$
Barometric pressure: 752.5 mmHg
Item functional (based on performance criteria):
All OK: Pass
Deviation from pretest: None

Third check:

Time: 1445
Temperature: $49^{\circ}\text{C} \pm 1^{\circ}\text{C}$
Humidity: $15\% \text{ RH} \pm 1\% \text{ RH}$
Barometric pressure: 752.5 mmHg
Item functional (based on performance criteria):
All OK: Pass
Deviation from pretest: None

POSTTEST DATA

Posttest performance check:

(complete check of item and accessories)

Time of test end: 1515

Item functional: (based on performance criteria):

All OK: Pass

Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: None

3.2.6 High Temperature Storage test

High Temperature test (Equipment in Storage) Report Form

Nomenclature: Molecular Sieve Oxygen Generating System
Manufacturer: Essex
Model number: 60C-0037-2
Serial number: 0002
Military item number: None

Options installed: None

Date of test: 25 Jun 92

Item configuration during test:

Sitting on wire test stand, not operating. The unit is in storage.

Performance test criteria:

The unit's ability to provide a continuous flow (16 lpm) of 90% (or higher) oxygen with an input pressure of 60 psi.

Ambient conditions outside chamber:

Temperature	28°C
Humidity	48% RH
Barometric pressure	752 mmHg

PRETEST DATA

Pretest performance check:

Item functional (based on performance criteria):
All OK: Pass

Installation of item in test facility:

list connections to power	None
list connections to simulators	None
list connections to dummy loads	None
list unconnected terminals	All
distance from north wall (meters)	0.56
distance from south wall (meters)	1.02
distance from east wall (meters)	1.57
distance from west wall (meters)	1.45
distance from ceiling (meters)	1.19
distance from floor (meters)	0.97

Time of test start: 0800
Mid-test time: 1100
Mid-test temperature: $71^{\circ}\text{C} \pm 1^{\circ}\text{C}$
Mid-test humidity: $15\% \text{ RH} \pm 1\% \text{ RH}$

POSTTEST DATA

Posttest performance check:
(complete check of item and accessories)
Time of test end: 1400
Item functional (based on performance criteria):
All OK: Pass
Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: None

3.2.7 Low Temperature Operation test

Low Temperature test (Equipment Operating) Report Form

Nomenclature: Molecular Sieve Oxygen Generating System
Manufacturer: Essex
Model number: 60C-0037-2
Serial number: 0002
Military item number: None

Options installed: None

Date of test: 22 Jun 92

Item configuration during test:

Operating on a wire test stand and connected to the
USAARL oxygen test apparatus via portals in the chamber
wall.

Performance test criteria:

The unit's ability to provide a continuous flow (16
lpm) of 90% (or higher) oxygen with an input pressure
of 60 psi.

Ambient conditions outside chamber:

Temperature	29°C
Humidity	66% RH
Barometric pressure	753 mmHg

PRETEST DATA

Pretest performance check:

Item functional (based on performance criteria)

All OK: Pass

Installation of item in test facility:

list connections to power	28 Vdc
list connections to simulators	Oxygen tester
list connections to dummy loads	None
list unconnected terminals	None

distance from north wall (meters)	0.56
distance from south wall (meters)	1.02
distance from east wall (meters)	1.57
distance from west wall (meters)	1.45
distance from ceiling (meters)	1.19
distance from floor (meters)	0.97

Time of test start: 0900

Performance checks during test:

First Check:

Time:	0930
Temperature:	- 25°C ± 1°C
Humidity:	NA
Barometric pressure:	753 mmHg
Item functional (based on performance criteria)	All OK: Pass
Deviation from pretest:	None

Second check:

Time:	1000
Temperature:	- 25°C ± 1°C
Humidity:	NA
Barometric pressure:	753 mmHg
Item functional (based on performance criteria)	All OK: Pass
Deviation from pretest:	None

Third check:

Time:	1030
Temperature:	- 25°C ± 1°C
Humidity:	NA
Barometric pressure:	753 mmHg
Item functional (based on performance criteria)	All OK: Pass
Deviation from pretest:	None

POSTTEST DATA

Posttest performance check:

(complete check of item and accessories)

Time of test end: 1100

Item functional (based on performance criteria)

All OK: Pass

Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: None

3.2.8 Low Temperature Storage test

Low Temperature test (Equipment in Storage) Report Form

Nomenclature: Molecular Sieve Oxygen Generating System
Manufacturer: Essex
Model number: 60C-0037-2
Serial number: 0002
Military item number: None

Options installed: None

Date of test: 26 Jun 92

Item configuration during test:

Sitting on wire test stand, not operating. The unit is in storage.

Performance test criteria:

The unit's ability to provide a continuous flow (16 lpm) of 90% (or higher) oxygen with an input pressure of 60 psi.

Ambient conditions outside chamber:

Temperature	29°C
Humidity	48% RH
Barometric pressure	753 mmHg

PRETEST DATA

Pretest performance check:

Item functional (based on performance criteria):
All OK: Pass

Installation of item in test facility:

list connections to power:	None
list connections to simulators:	None
list connections to dummy loads	None
list unconnected terminals	All
distance from north wall (meters)	0.56
distance from south wall (meters)	1.02
distance from east wall (meters)	1.57
distance from west wall (meters)	1.45
distance from ceiling (meters)	1.19
distance from floor (meters)	0.97

Time of test start: 0800
Mid-test time: 1100
Mid-test temperature: $-46^{\circ}\text{C} \pm 1^{\circ}\text{C}$

POSTTEST DATA

Posttest performance check:
(complete check of item and accessories)
Time of test end: 1400
Item functional (based on performance criteria):
All OK: Pass
Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: None

3.2.9 Humidity test

Humidity Test Report Form

Nomenclature: Molecular Sieve Oxygen Generating System
Manufacturer: Essex
Model number: 60C-0037-2
Serial number: 0002
Military item number: None

Options installed: None

Date of test: 24 Jun 92

Item configuration during test:

Operating on a wire test stand and connected to the
USAARL oxygen test apparatus via portals in the chamber
wall.

Performance test criteria:

The unit's ability to provide a continuous flow (27
lpm) of 90% (or higher) oxygen with an input pressure
of 60 psi.

Ambient conditions outside chamber:

Temperature	29°C
Humidity	48% RH
Barometric pressure	754 mmHg

PRETEST DATA

Pretest performance check:

Item functional (based on performance criteria):
all OK: Pass

Installation of item in test facility:

list connections to power	28 Vdc
list connections to simulators	OTA
list connections to dummy loads	None
list unconnected terminals	None
distance from north wall (meters)	0.56
distance from south wall (meters)	1.02
distance from east wall (meters)	1.57
distance from west wall (meters)	1.45
distance from ceiling (meters)	1.19
distance from floor (meters)	0.97

IN-TEST DATA

Time of test start: 0900

Performance checks during test:

First check:

Time: 0945
Temperature: $29.5^{\circ}\text{C} \pm 1^{\circ}\text{C}$
Humidity: $95\% \text{ RH} \pm 1\% \text{ RH}$
Barometric pressure: 754 mmHg
Item functional (based on performance criteria):
All OK: Pass
Deviation from pretest: none

Second check:

Time: 1030
Temperature: $29.5^{\circ}\text{C} \pm 1^{\circ}\text{C}$
Humidity: $95\% \text{ RH} \pm 1\% \text{ RH}$
Barometric pressure: 754 mmHg
Item functional (based on performance criteria):
All OK: Pass
Deviation from pretest: None

Third check:

Time: 1115
Temperature: $29.5^{\circ}\text{C} \pm 1^{\circ}\text{C}$
Humidity: $95\% \text{ RH} \pm 1\% \text{ RH}$
Barometric pressure: 754 mmHg
Item functional (based on performance criteria):
All OK: Pass
Deviation from pretest: None

Fourth check:

Time: 1200
Temperature: $29.5^{\circ}\text{C} \pm 1^{\circ}\text{C}$
Humidity: $95\% \text{ RH} \pm 1\% \text{ RH}$
Barometric pressure: 754 mmHg
Item functional (based on performance criteria):
All OK: Pass
Deviation from pretest: None

Fifth check:

Time: 1245
Temperature: 29.5°C \pm 1°C
Humidity: 95% RH \pm 1% RH
Barometric pressure: 754 mmHg
Item functional (based on performance criteria):
All OK: Pass
Deviation from pretest: None

POSTTEST DATA

Posttest performance check:
(complete check of item and accessories)
Time of test end: 1300
Item functional (based on performance criteria):
All OK: Pass
Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: None

3.2.10

Electromagnetic Characteristics testing
Evaluation of Performance

Nomenclature: Molecular Sieve Oxygen Generating System
Manufacturer: Essex
Model number: 60C-0037-2
Serial number: 0002
Military item number: none

Options installed: None

Date of test: 3 Aug 92

<u>Tests</u>	<u>Testing configuration(s):</u>	<u>Performance</u> <u>(Pass/fail)</u>
<u>Conducted Emissions</u>	Operating on copper bench	Fail
CE01		
Comments:	Failure emissions detected in the following range. 23.42 - 50 MHz with 0.1 - 7.9 dB of failure	
CE02		NA
Comments:	Unit does not have ac power leads.	
CE04	Operating on copper bench.	Fail
Comments:	Broadband emission failures detected in the following range: 0.325 - 45 MHz with 0.1 - 44.9 dB of failure.	
<u>Conducted Susceptibility</u>		
CS02	Operating on copper bench.	Not conducted
Comments:	Not able to measure test generated signals on power conductors.	
CS06	Operating on copper test bench.	Pass
Comments:	Not susceptible to test spikes.	

Radiated
Emissions

RE02 Operating on wooden test stand. Fail
Comments: Narrowband and broadband emissions were detected.
 NB: 3.789 MHz with 52.7 dB of failure
 BB: .014 - 625 MHz with 0.4 dB - 32.6 dB of
 failure.

Radiated
Susceptibility

RS03 Operating on wooden test stand. Pass
Comments: Not susceptible to test generated signals.

3.2.11 Acoustical test

Acoustical test Report Form

Nomenclature: Molecular Sieve Oxygen Generating System
 Manufacturer: Essex
 Model number: 60C-0037-2
 Serial number: 0002
 Military item number: Ncne

Options installed: None

Date of test: 23 Jun 92

1/3 Octave
 Band Center
 Frequency
 (Hz)

Sound Pressure Level
 (dB)

	Orientation			
	East	North	South	West
20	65.6	69.1	61.8	61.8
30	67	68.9	63	61.8
40	67.6	67.5	63.9	62.6
50	67.8	67.3	65.5	64.9
60	68.6	69.6	70.5	66.7
70	67.8	67.5	65.3	65.7
80	67.4	67.1	68.6	67.9
90	66.3	65.9	65	64.4
100	62.7	63.6	61	61.9
200	62	61.8	58.9	59
300	60.6	61.4	58.3	60
400	56.5	57.2	55.2	55.1
500	56.7	58.5	54.1	57
600	60.1	63.1	58.6	61.5
700	71	73.3	72.5	72.6
800	72.8	72.8	74.3	73.4
900	60.6	64.8	61.2	59.2
1000	71.4	69.9	70.6	68
1100	73.9	69.9	73.7	70.4
1200	70.7	73.2	66.9	70.3
1300	75.7	80.6	68	72.7
1400	75.7	76.2	67	71.2
1500	73.7	70.2	68.8	70.2
1600	77.6	75.7	68.8	73.2
1700	73.3	69.7	64.6	68.5
1800	75	70.4	66.2	69.4
1900	74.7	66.1	66.1	70
10000	73.7	69.5	65.1	66.3

1/3 Octave
Band Center
Frequency
(Hz)

Sound Pressure Level
(dB)

Orientation

East

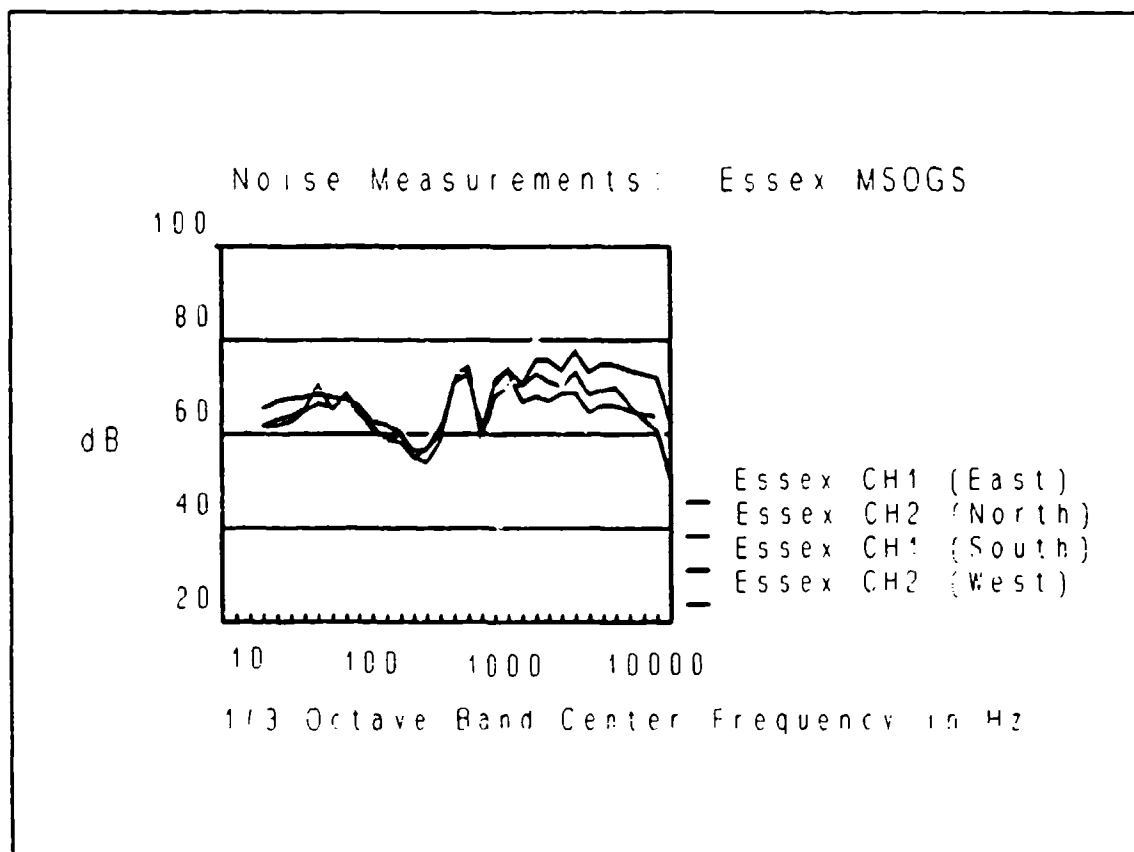
North

South

West

11000	72.7	69	64.1	63.1
12000	72	63.9	63.7	60.8
13000	62.3	54.2	53.7	50.5

Note: A graphical representation of test results follows:



3.3 CRITERIA, SIGNIFICANT PROBLEMS, AND SUGGESTED IMPROVEMENTS

3.3.1 Criteria

<u>Item no.</u>	<u>Criteria (source)</u>	<u>Remarks</u>	<u>Applicable subparagraph</u>
1	The physical inventory is conducted solely for investigation and documentation.	NA	2.1.2.1
2	The Essex Model 60C-0037-2 will perform consistently and accurately with an output flow of 16 lpm.	met	2.1.2.2
3	The Essex Model 60C-0037-2 shall meet the standards established by NAEP 99 for electrical safety of medical equipment.	partially met	2.2.2
4	The unit will produce oxygen of at least 90% concentration at the manufacturer's specified rate.	met	2.3.2
5	The MSOGS unit must be rated satisfactory in all major categories of the human factors evaluation.	met	2.4.2
6	The Essex Model 60C-0037-2 will demonstrate proper operation - while exposed to an altitude equivalency of 15,000 feet above sea level.	met	2.5.2
7	The Essex Model 60C-0037-2 will remain operational while exposed to vibrational stresses.	met	2.6.2
8	The Essex Model 60C-0037-2 will remain operational during the high temperature operation check.	met	2.7.2.1
9	The Essex Model 60C-0037-2 will remain operational after the high temperature storage.	met	2.7.2.2

10	The Essex Model 60C-0037-2 will remain operational during the low temperature operation check.	met	2.8.2.1
11	The Essex Model 60C-0037-2 will remain operational after the low temperature storage.	met	2.8.2.2
12	The Essex Model 60C-0037-2 will remain operational while exposed to a high humidity.	met	2.9.2
13	The Essex Model 60C-0037-2 will not produce emissions in excess of the limits set forth in MIL-STD-461C, Notice 2.	par- tially met	2.10.2.1
14	The Essex Model 60C-0037-2 will not malfunction when it is subjected to radiated fields as specified in MIL-STD-461C, Notice 2.	met	2.10.2.2
15	The Essex Model 60C-0037-2 will not conduct emissions in excess of the limits set forth in MIL-STD-461C, Notice 2.	par- tially met	2.10.2.3
16	The Essex Model 60C-0037-2 will not malfunction when it is subjected to conducted emissions as specified in MIL-STD-461C, Notice 2.	met	2.10.2.4
17	The Essex MSOGS will not produce a noise hazard during normal operation.	not evalu- ated	2.11.2

3.3.2 Significant problems which require corrective action

None

3.3.3 Suggested improvements

None

3.4 REFERENCES

3.4.1 Department of Defense. 1986. EMI characteristics, requirements for equipment. Washington, DC. MIL-STD-461C, Notice 2. August.

3.4.2 Department of Defense. 1971. EMI characteristics, measurement of. Washington, DC. MIL-STD-462, Notice 5. February.

3.4.3 Department of Defense. 1983. Environmental test methods and engineering guidelines. Washington, DC. MIL-STD-810E. July.

3.4.4 Department of the Army. 1987. Maintenance management procedures for medical equipment. Washington, DC. TB 38-750-2. April.

3.4.5 National Fire Protection Association. 1987. Standard for health care facilities. Quincy, MA. NAF 99. February.

3.5 ABBREVIATIONS

ac	alternating current
AVSCOM	U.S. Army Aviation Systems Command
AEST	aeromedical equipment suitability test
AWR	airworthiness release
BB	broadband
BPM	beats per minute
CAAF	Cairns Army Airfield
CRT	cathode ray tube
dB	decibel
dc	direct current
ECG	electrocardiograph
EMC	electromagnetic compatibility
EMI	electromagnetic interference
fpm	feet per minute
GFE	government furnished equipment
Gpk	gravity, peak
G(rms)	gravity (root mean square)
Hz	hertz
IAW	in accordance with
ITOP	in-flight test operating procedure
IGE	in-ground effect
kHz	kilohertz
KIAS	knots indicated airspeed
LCD	liquid crystal display
LED	light emitting diode
LISN	line impedance stabilization networks
MEDEVAC	medical evacuation
MHz	mega hertz
MIL-STD	military standard
mL	milliliter
mm	millimeter
mmHg	millimeters of mercury
m.s.l.	mean sea level
NAFP	National Association of Fire Prevention
NB	narrowband
NBC	nuclear, biological, and chemical
NiCad	nickel cadmium

NVG	night vision goggles
oct	octave
RAM	random access memory
RF	radio frequency
RH	relative humidity
ROM	read only memory
TB	technical bulletin
TFT	technical feasibility testing
T & E	test and evaluation
UES	Universal Energy Systems, Inc.
USAARL	U.S. Army Aeromedical Research Laboratory
V/m	volts per meter

3.6 LIST OF MANUFACTURERS

- 3.6.1 Essex Instruments and Life Support
P.O. Box 6234
Fort Worth, TX 76115
- 3.6.2 Neurodyne-Dempsey, Inc.
200 Arrowhead Drive
Carson City, NV 89701
- 3.6.3 Tenney Engineering, Inc.
1090 Springfield Road
Post Office Box 3142
Union, NJ 07083
- 3.6.4 Unholtz-Dickey Corporation
6 Brookside Drive
Wallingford, CT 06492
- 3.6.5 Solar Electronics Company
901 North Highland Avenue
Hollywood, CA 90038
- 3.6.6 Tektronix, Inc.
P.O. Box 500
Beaverton, OR 97077
- 3.6.7 Nagra Magnetics, Inc.
19 West 44th Street
Room 715
New York, NY 10036

3.7 INITIAL DISTRIBUTION

Commander, U.S. Army Natick Research,
Development and Engineering Center
ATTN: SATNC-MIL (Documents
Librarian)
Natick, MA 01760-5040

U.S. Army Communications-Electronics
Command
ATTN: AMSEL-RD-ESA-D
Fort Monmouth, NJ 07703

Commander
10th Medical Laboratory
ATTN: Audiologist
APO New York 09180

Naval Air Development Center
Technical Information Division
Technical Support Detachment
Warminster, PA 18974

Commanding Officer, Naval Medical
Research and Development Command
National Naval Medical Center
Bethesda, MD 20814-5044

Deputy Director, Defense Research
and Engineering
ATTN: Military Assistant
for Medical and Life Sciences
Washington, DC 20301-3080

Commander, U.S. Army Research
Institute of Environmental Medicine
Natick, MA 01760

Library
Naval Submarine Medical Research Lab
Box 900, Naval Sub Base
Groton, CT 06349-5900

Director, U.S. Army Human
Engineering Laboratory
ATTN: Technical Library
Aberdeen Proving Ground, MD 21005

Commander
Man-Machine Integration System
Code 602
Naval Air Development Center
Warminster, PA 18974

Commander
Naval Air Development Center
ATTN: Code 602-B (Mr. Brindle)
Warminster, PA 18974

Commanding Officer
Armstrong Laboratory
Wright-Patterson
Air Force Base, OH 45433-6573

Director
Army Audiology and Speech Center
Walter Reed Army Medical Center
Washington, DC 20307-5001

Commander/Director
U.S. Army Combat Surveillance
and Target Acquisition Lab
ATTN: DELCS-D
Fort Monmouth, NJ 07703-5304

Commander, U.S. Army Institute
of Dental Research
ATTN: Jean A. Setterstrom, Ph. D.
Walter Reed Army Medical Center
Washington, DC 20307-5300

Commander, U.S. Army Test
and Evaluation Command
ATTN: AMSTE-AD-H
Aberdeen Proving Ground, MD 21005

Naval Air Systems Command
Technical Air Library 950D
Room 278, Jefferson Plaza II
Department of the Navy
Washington, DC 20361

Director
U.S. Army Ballistic
Research Laboratory
ATTN: DRXBR-OD-ST Tech Reports
Aberdeen Proving Ground, MD 21005

Commander
U.S. Army Medical Research
Institute of Chemical Defense
ATTN: SGRD-UV-AO
Aberdeen Proving Ground,
MD 21010-5425

Commander, U.S. Army Medical
Research and Development Command
ATTN: SGRD-RMS (Ms. Madigan)
Fort Detrick, Frederick, MD 21702-5012

Director
Walter Reed Army Institute of Research
Washington, DC 20307-5100

HQ DA (DASG-PSP-O)
5109 Leesburg Pike
Falls Church, VA 22041-3258

Harry Diamond Laboratories
ATTN: Technical Information Branch
2800 Powder Mill Road
Adelphi, MD 20783-1197

U.S. Army Materiel Systems
Analysis Agency
ATTN: AMXSY-PA (Reports Processing)
Aberdeen Proving Ground
MD 21005-5071

U.S. Army Ordnance Center
and School Library
Simpson Hall, Building 3071
Aberdeen Proving Ground, MD 21005

U.S. Army Environmental
Hygiene Agency
Building E2100
Aberdeen Proving Ground, MD 21010

Technical Library Chemical Research
and Development Center
Aberdeen Proving Ground, MD
21010-5423

Commander
U.S. Army Medical Research
Institute of Infectious Disease
SGRD-UIZ-C
Fort Detrick, Frederick, MD 21702

Director, Biological
Sciences Division
Office of Naval Research
600 North Quincy Street
Arlington, VA 22217

Commander
U.S. Army Materiel Command
ATTN: AMCDE-XS
5001 Eisenhower Avenue
Alexandria, VA 22333

Commandant
U.S. Army Aviation
Logistics School ATTN: ATSQ-TDN
Fort Eustis, VA 23604

Headquarters (ATMD)
U.S. Army Training
and Doctrine Command
ATTN: ATBO-M
Fort Monroe, VA 23651

IAF Liaison Officer for Safety
USAF Safety Agency/SEFF
9750 Avenue G, SE
Kirtland Air Force Base
NM 87117-5671

Naval Aerospace Medical
Institute Library
Building 1953, Code 03L
Pensacola, FL 32508-5600

Command Surgeon
HQ USCENTCOM (CCSG)
U.S. Central Command
MacDill Air Force Base, FL 33608

Air University Library
(AUL/LSE)
Maxwell Air Force Base, AL 36112

U.S. Air Force Institute
of Technology (AFIT/LDEE)
Building 640, Area B
Wright-Patterson
Air Force Base, OH 45433

Henry L. Taylor
Director, Institute of Aviation
University of Illinois-Willard Airport
Savoy, IL 61874

Chief, National Guard Bureau
ATTN: NGB-ARS (COL Urbauer)
Room 410, Park Center 4
4501 Ford Avenue
Alexandria, VA 22302-1451

Commander
U.S. Army Aviation Systems Command
ATTN: SGRD-UAX-AL
4300 Goodfellow Blvd., Building 105
St. Louis, MO 63120

U.S. Army Aviation Systems Command
Library and Information Center Branch
ATTN: AMSAV-DIL
4300 Goodfellow Boulevard
St. Louis, MO 63120

Federal Aviation Administration
Civil Aeromedical Institute
Library AAM-400A
P.O. Box 25082
Oklahoma City, OK 73125

Commander
U.S. Army Academy
of Health Sciences
ATTN: Library
Fort Sam Houston, TX 78234

Commander
U.S. Army Institute of Surgical Research
ATTN: SGRD-USM (Jan Duke)
Fort Sam Houston, TX 78234-6200

AAMRL/HEX
Wright-Patterson
Air Force Base, OH 45433

John A. Dellinger,
Southwest Research Institute
P. O. Box 28510
San Antonio, TX 78284

Product Manager
Aviation Life Support Equipment
ATTN: AMCPM-ALSE
4300 Goodfellow Boulevard
St. Louis, MO 63120-1798

Commander
U.S. Army Aviation
Systems Command
ATTN: AMSAV-ED
4300 Goodfellow Boulevard
St. Louis, MO 63120

Commanding Officer
Naval Biodynamics Laboratory
P.O. Box 24907
New Orleans, LA 70189-0407

Assistant Commandant
U.S. Army Field Artillery School
ATTN: Morris Swott Technical Library
Fort Sill, OK 73503-0312

Mr. Peter Seib
Human Engineering Crew Station
Box 266
Westland Helicopters Limited
Yeovil, Somerset BA20 2YB UK

U.S. Army Dugway Proving Ground
Technical Library, Building 5330
Dugway, UT 84022

U.S. Army Yuma Proving Ground
Technical Library
Yuma, AZ 85364

AFFTC Technical Library
6510 TW/TSTL
Edwards Air Force Base,
CA 93523-5000

Commander
Code 3431
Naval Weapons Center
China Lake, CA 93555

Aeromechanics Laboratory
U.S. Army Research and Technical Labs
Ames Research Center, M/S 215-1
Moffett Field, CA 94035

Sixth U.S. Army
ATTN: SMA
Presidio of San Francisco, CA 94129

Commander
U.S. Army Aeromedical Center
Fort Rucker, AL 36362

Strughold Aeromedical Library
Document Service Section
2511 Kennedy Circle
Brooks Air Force Base, TX 78235-5122

Dr. Diane Damos
Department of Human Factors
ISSM, USC
Los Angeles, CA 90089-0021

U.S. Army White Sands
Missile Range
ATTN: STEWS-IM-ST
White Sands Missile Range, NM 88002

U.S. Army Aviation Engineering
Flight Activity
ATTN: SAVTE-M (Tech Lib) Stop 217
Edwards Air Force Base, CA 93523-5000

Ms. Sandra G. Hart
Ames Research Center
MS 262-3
Moffett Field, CA 94035

Commander, Letterman Army Institute
of Research
ATTN: Medical Research Library
Presidio of San Francisco, CA 94129

Commander
U.S. Army Medical Materiel
Development Activity
Fort Detrick, Frederick, MD 21702-5009

Commander
U.S. Army Health Services Command
ATTN: HSOP-SO
Fort Sam Houston, TX 78234-6000

U. S. Army Research Institute
Aviation R&D Activity
ATTN: PERI-IR
Fort Rucker, AL 36362

Commander
U.S. Army Safety Center
Fort Rucker, AL 36362

U.S. Army Aircraft Development
Test Activity
ATTN: STEBG-MP-P
Cairns Army Air Field
Fort Rucker, AL 36362

Commander, U.S. Army Medical Research
and Development Command
ATTN: SGRD-PLC (COL Schnakenberg)
Fort Detrick, Frederick, MD 21702

TRADOC Aviation LO
Unit 21551, Box A-209-A
APO AE 09777

Netherlands Army Liaison Office
Building 602
Fort Rucker, AL 36362

British Army Liaison Office
Building 602
Fort Rucker, AL 36362

Italian Army Liaison Office
Building 602
Fort Rucker, AL 36362

Directorate of Training Development
Building 502
Fort Rucker, AL 36362

Chief
USAHEL/USAAVNC Field Office
P. O. Box 716
Fort Rucker, AL 36362-5349

Commander, U.S. Army Aviation Center
and Fort Rucker
ATTN: ATZQ-CG
Fort Rucker, AL 36362

Chief
Test & Evaluation Coordinating Board
Cairns Army Air Field
Fort Rucker, AL 36362

MAJ Terry Newman
Canadian Army Liaison Office
Building 602
Fort Rucker, AL 36362

German Army Liaison Office
Building 602
Fort Rucker, AL 36362

LTC Patrice Cottebrune
French Army Liaison Office
USAAVNC (Building 602)
Fort Rucker, AL 36362-5021

Australian Army Liaison Office
Building 602
Fort Rucker, AL 36362

Dr. Garrison Rapmund
6 Burning Tree Court
Bethesda, MD 20817

Commandant, Royal Air Force
Institute of Aviation Medicine
Farnborough, Hampshire GU14 6SZ UK

Commander
U.S. Army Biomedical Research
and Development Laboratory
ATTN: SGRD-UBZ-I
Fort Detrick, Frederick, MD 21702

Defense Technical Information
Cameron Station, Building 5
Alexandra, VA 22304-6145

Commander, U.S. Army Foreign Science
and Technology Center
AIFRTA (Davis)
220 7th Street, NE
Charlottesville, VA 22901-5396

Director,
Applied Technology Laboratory
USARTL-AVSCOM
ATTN: Library, Building 401
Fort Eustis, VA 23604

Commander, U.S. Air Force
Development Test Center
101 West D Avenue, Suite 117
Eglin Air Force Base, FL 32542-5495

Aviation Medicine Clinic
TMC #22, SAAF
Fort Bragg, NC 28305

Dr. H. Dix Christensen
Bio-Medical Science Building, Room 753
Post Office Box 26901
Oklahoma City, OK 73190

Commander, U.S. Army Missile
Command
Redstone Scientific Information Center
ATTN: AMSMI-RD-CS-R
/ILL Documents
Redstone Arsenal, AL 35898

Director
Army Personnel Research Establishment
Farnborough, Hants GU14 6SZ UK

U.S. Army Research and Technology
Laboratories (AVSCOM)
Propulsion Laboratory MS 302-2
NASA Lewis Research Center
Cleveland, OH 44135

COL John F. Glenn
U.S. Army Medical Research
& Development Command
SGRD-ZC
Fort Detrick, Frederick, MD 21702-5012

Dr. Eugene S. Channing
7985 Schooner Court
Frederick, MD 21701-3273

USAMRDC Liaison at Academy
of Health Sciences
ATTN: HSHA-ZAC-F
Fort Sam Houston, TX 78234

Dr. A. Kornfield, President
Biosearch Company
3016 Revere Road
Drexel Hill, PA 29026

NVEOD
AMSEL-RD-ASID
(Attn: Trang Bui)
Fort Belvoir, VA 22060

CA Av Med
HQ DAAC
Middle Wallop
Stockbridge, Hants S020 8DY UK

Commander and Director
USAE Waterways Experiment Station
ATTN: CEWES-IM-MI-R
Alfrieda S. Clark, CD Dept.
3909 Halls Ferry Road
Vicksburg, MS 39180-6199

Dr. Christine Schlichting
Behavioral Sciences Department
Box 900, NAVUBASE NLON
Groton, CT 06349-5900

Brazilian Army Liaison Office
Building 602
Fort Rucker, AL 36362

COL C. Fred Tyner
U.S. Army Medical Research
& Development Command
SGRD-ZB
Fort Detrick, Frederick, MD 21702-5012

Director
Directorate of Combat Developments
ATZQ-CD
Building 515
Fort Rucker, AL 36362

**END
FILMED**

DATE:

10-93

DTIC